

Sterile SatPax® MicroSeal® WFI

Sterile Pre-wetted Sealed Edge Cleanroom Laundered 100% Polyester Knit Wiper

Sterile SatPax® MicroSeal® WFI combines MicroSeal® 1200 ultrasonically sealed edge cleanroom laundered wipers composed of 100% continuous filament polyester knit fabric with a high saturation level of 70% IPA and 30% WFI water.

This pre-wetted format provides a cost effective and easy to use solution versus traditional bulk handling of solvents, maintenance of squirt bottles and inconsistent wetting and cleaning associated with wetting a dry wiper.

Each package is then gamma irradiated sterile to a 10^{-6} Sterility Assurance Level making the product ideal for cleaning critical sterile processing environments.



Key Attributes

- 100% continuous filament polyester knit
- Ultrasonically sealed edge for reduced fiber contamination
- Laundered and packaged in Berkshire's ISO Class 4 cleanroom
- Pre-wetted with consistent 70% IPA/ 30% WFI (Water for Injection) to a high saturation level
- Re-sealable solvent resistant packaging
- Triple packaged for easy entry into sterile environments
- Gamma irradiated and sterile validated to a 10^{-6} Sterility Assurance Level per AAMI guidelines
- Certificate of Sterility with each case and available on the web by lot number 24/7. Details: Expiration Date, Radiation Dosage, Lot Information, Meets cGMP requirements for traceability

Benefits

- Critically low fibers, particles, ions and extractables combined with high abrasion resistance
- WFI (Water for Injection) solution for critical applications requiring WFI
- Reduces alcohol usage and preparation / handling costs
- Reduces VOC emissions
- Increases cleaning efficiency
- Increases cleaning protocol consistency

Environmental

- Reduces VOC (Volatile Organic Compound) emissions
- Complies with RoHS and pre-registered under REACH

Applications

- Designed for use in ISO Class 3 and higher sterile cleanroom environments
- Designed for use in wet cleaning of critical surfaces where control of flammable solvents and flammable solvent concentrations is required
- Final cleaning of surfaces or products prior to manufacturing or packaging
- High saturation level is ideal for removing cleaning and disinfecting residues in sterile environments

Saturation Levels

The amount of solution contained in each wiper can be varied according to customer requirements. Higher saturation levels apply more solution to the surface during cleaning.

Validated Sterile

ANSI/AAMI/ISO procedures have been followed for verification of irradiation dose with audits being performed to monitor bioburden levels and dose. The Sterility Assurance Level (SAL) or the probability of a viable organism being present on a product unit after sterilization is a SAL of 10^{-6} (which is the chance of 1 non-sterile unit in 1,000,000).

Other Sterile Pre-wetted wipers

- Sterile SatPax® 670
- Sterile SatPax® 1000
- Sterile SatPax® 1200

www.berkshire.com.sg

Contact: Tel 65.6252.4313

Fax: 65.6252.4312

Email:enquiries@berkshire.com.sg

Americas	Tel 1 413 528 2602	Fax 1 413 528 2614	info@berkshire.com
Europe	Tel 44 (0) 870 757 2877	Fax 44 (0) 870 757 2878	enquiries@berkshire.uk.com
SE Asia	Tel 65 6252 4313	Fax 65 6252 4312	enquiries@berkshire.com.sg
Japan	Tel 81 3 5827 2380	Fax 81 3 5827 2382	master@berkshire.co.jp

Technical Data (In Dry State)

Attribute		Units	Value	Test Method
Basis Weight		g/m ²	154	TAPPI T-410
Caliper		µm	507	TAPPI T-411
Fibers	≥100µm	fibers/cm ²	0.036	TEST-RP.CC004.3, Sec 6.1.3 / Sec 6.2.2
Particles	≥0.5µm	x10 ³ /cm ²	0.40	TEST-RP.CC004.3, Sec 6.1.3 / Sec 6.2.1
Sorbency	Capacity	mL/m ²	512	TEST-RP.CC004.3, Sec 8.1 modified / Sec 8.2 modified
	Efficiency	mL/g	3.3	
	Rate	seconds	1	
Non-Volatile Residue	DI Water	g/m ²	0.045	TEST-RP.CC004.3, Sec 7.1.2
	IPA	g/m ²	0.0082	
Ions	Na ⁺	ppm	0.17	TEST-RP.CC004.3, Sec 7.2.2
	K ⁺	ppm	0.010	
	Ca ⁺⁺	ppm	0.046	
	Mg ⁺⁺	ppm	0.011	
	Cl ⁻	ppm	0.30	

Notes:

- Technical data represented in this table are typical values at the time of publication. These should not be used as product specifications.
- Due to differences in test methods applied and equipment utilized by different wiper manufacturers, valid product comparisons may only be obtained through side-by-side testing in the same test facility, under similar conditions
- Third party testing can be performed upon request

Order Information:

Product	Number	Size	Shts/pk	Pks/cs	IPA/WFI	Saturation	VOC % by	Style
Sterile SatPax® MS-WFI	SSPMS015WFI	12x12" (30x30cm)	50	5	70/30	60%	42%	Flat Pack

Other Berkshire Products



Wipers




Mop Systems



Face Masks



Glove Liners

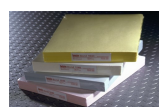

Documentation Systems



Swabs

www.berkshire.com.sg

Contact: Tel 65.6252.4313 / Fax: 65.6252.4312

Email: enquiries@berkshire.com.sg

BC-3172_08-2015